

WHY SPRITZTUBE® SHOULD BE USED

OXYGENATE QUICKLY	WORKS FOR LOW-SKILL OPERATORS	BOTH FOR EMERGENCY AND HOSPITAL SETTINGS	CAN BE INSERTED BLINDLY
GIVES TIME TO THINK ABOUT FURTHER STEPS WHILE VENTILATING THE PATIENT	WORKS AS ENDOTRACHEAL TUBE DURING INTUBATION – NO NEED FOR EXTRA MANEUVERES TO EXCHANGE AIRWAY DEVICES	INTUBATION IS POSSIBLE WITH ANY EQUIPMENT AVAILABLE	IDEAL COMBINED WITH FLEXIBLE FIBERSCOPE INTUBATION – APNEA FASE REDUCED TO MINIMUM

TECHNICAL DATASHEET


Materials: Tube - medical grade silicone with embedded metal coil
 Balloons - medical grade silicone
 Cuffs - medical grade silicone
 Stylet - silicone coated malleable steel

Usage: Resterilizable in autoclave up to 30 times (tracked by a progressive code number on the tube)

Safety: Latex free

Certification: CE marked (93/42/EEC Directive)
 Certificate MED 9827-A by KIWA CERMET ITALIA spa (Notified body N.0476)

Classification: Class IIa - Rule 5

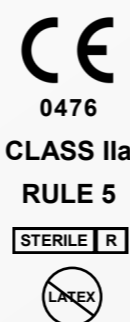

ORDERING INFORMATION

REF	CH	ext. Ø	int. Ø	Shore	Length	Packaging	Box	Shelf life
MED ST-32	32	10,6 mm	7,5 mm	50	330 mm	Single sterile pouch, sterilized by gamma rays	3 pcs	5 years
MED ST-28	28	9,3 mm	6,5 mm	50	330 mm	Single sterile pouch, sterilized by gamma rays	3 pcs	5 years



Manufactured by:



Distributed by:



Via Gramsci, 9 - 40018 San Pietro in Casale (BO) - Italy
 Tel. +39 051.817902 +39 051.817850 Fax +39 051.817845
 www.simitalia.it www.medeurope.eu
 E-Mail: info@medeurope.eu


SPRITZTUBE®

by MED EUROPE s.r.l.



NEW TRACHEAL TUBE SIMPLIFYING THE ANESTHESIOLOGIST'S JOB

Patented device for **advanced and complete airway management**, that consists of a unique tube in medical grade silicone, featuring a malleable stylet with non-traumatic tip and **two separate cuffs, one supraglottic and one esophageal**, each having a dedicated inflation line

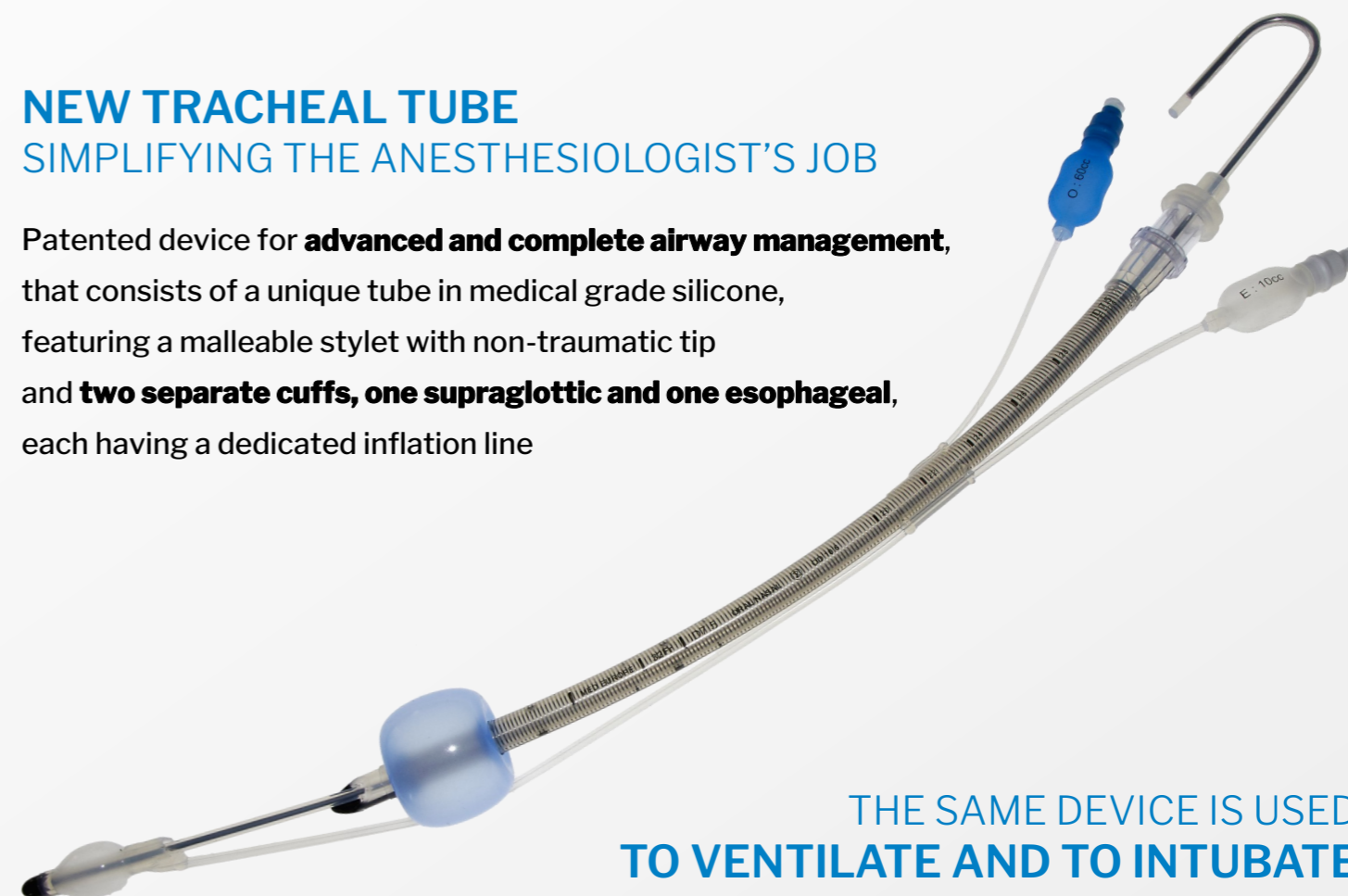

THE SAME DEVICE IS USED
TO VENTILATE AND TO INTUBATE

The Spritztube® allows rapid oxygenation and supraglottic ventilation, as well as subsequent **intubation with the aid of any available device (classic or video laryngoscope, rigid or flexible fiberscope)**, minimizing the apnea phase

EASY AND SAFE TO USE EVEN IN OUT-OF-HOSPITAL SETTING


Once inserted blindly or using a laryngoscope, even by nursing staff, the Spritztube® ensures the patients' ventilation until their arrival in hospital setting. Then, if intubation is required, it allows to proceed with it directly, simply by modifying its configuration.

No need for additional maneuvers to exchange the device!



WHAT IS SPRITZTUBE®

Magic in simplicity: **independent cuffs!**

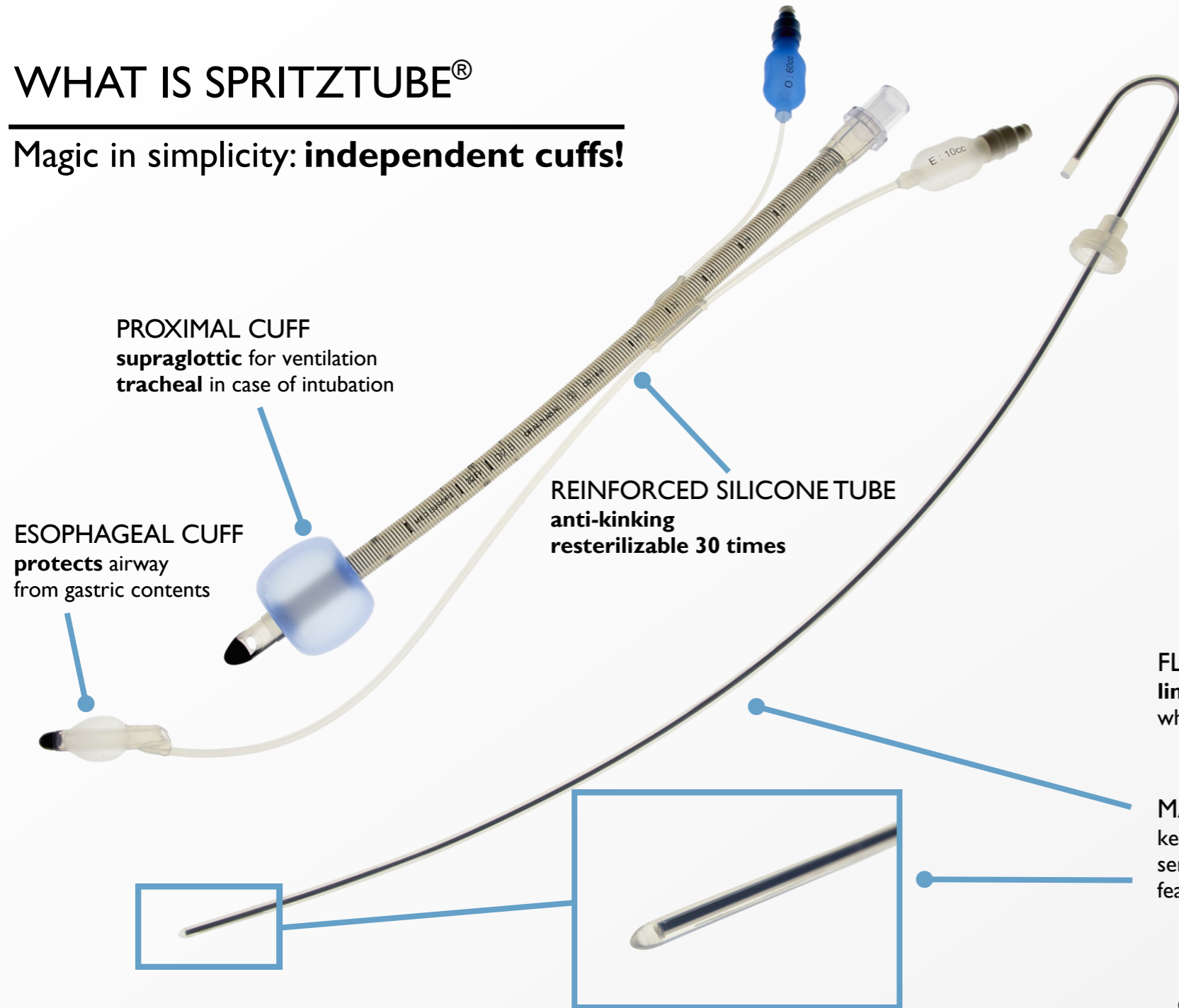

Spritztube® is an extraglottic airway device designed to ventilate as well as to permit tracheal intubation, increasing ease of use and success rate for the operators and improving safety of the patients, both in emergency settings and in elective surgery.

It's possible thanks to the patented flexibility of its design. It consists of a unique silicone tube with two separate cuffs, one proximal (oropharyngeal) and one distal (esophageal), linked to the tube by a flexible tab and easily removable when not needed.

With both cuffs inflated, Spritztube® works as an extraglottic airway device and ensures the patient's ventilation, giving the anaesthesiologist time to evaluate the next steps.

If the intubation is required, the proximal cuff is deflated, and Spritztube® is advanced through the glottis with the laryngoscope guidance or thanks to the use of fiberscope, which can be inserted directly through the Spritztube®, without stopping the ventilation.

Any existing intubation technique is possible, depending on the operator's preference and skill level. This implicates that the intubation maneuvers are simplified, since they don't require additional external endotracheal tube or tube exchange procedures.

In fact Spritztube® is the only existing device capable of turning from extraglottic to intraglottic because it itself acts as an endotracheal tube.

FLEXIBLE TAB
links the two cuffs while keeping them **independent**

MALLEABLE STYLET
keeps the cuffs in line during introduction serves as a guide during intubation features a **non traumatic tip**

SPRITZTUBE® is the answer
VENTILATE & INTUBATE with the same device

AIRWAY MANAGEMENT CHALLENGES

