


ByBand

Ring for bypass / sleeve / one anastomosis gastric bypass

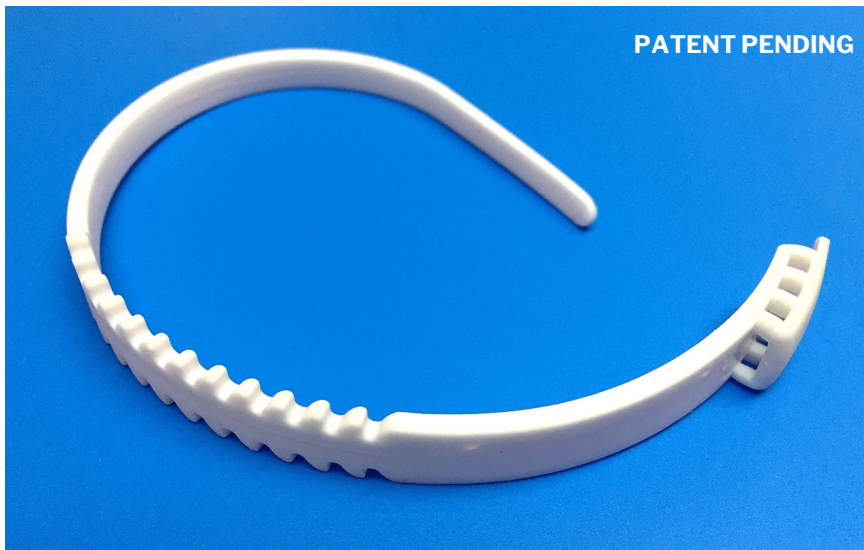

The device is intended to be used in the procedures of bariatric surgery. It can be placed around the pouch in gastric bypass, or at the upper part of sleeve gastrectomy and in mini-gastric bypass, with the intent to prevent their dilation. It is indicated for both primary and revisional bariatric surgery.

The weight regain over time, in bypass or sleeve gastrectomy is often linked to a gradual gastric pouch dilation. The purpose of the ByBand placement is to maintain the restrictive component in the gastric bypass and in the sleeve gastrectomy, over time.

The device is a pre-shaped implantable silicone ring. It is designed to have different locking positions that correspond to different diameters. ByBand can be implanted with either laparoscopic or laparotomic technique.

Characteristics

- Section measures: 13 x 3 mm.
- Length: 270 mm.
- Circumference corresponding to internal diameter: see the table.
- Material: medical grade implantable silicone, white color, x-ray, shore 70 A.
- Sterilised by gamma rays.
- Latex free.

CIRCUMF. mm	INTERNAL DIAM. mm
65	20,7
70	22,3
75	23,9
80	25,5
85	27
90	28,6
95	30,2
100	31,8
105	33,4
110	35
115	36,6
120	38,2

REF BB 05-07 L

Packaging: 1 piece/box

Manufactured by:



Via Gramsci, 9 - 40018 San Pietro in Casale (BO) - Italy
Tel. +39 051.817902 +39 051.817850 Fax +39 051.817845
www.simitalia.it www.medeurope.eu
E-Mail: info@medeurope.eu

Distributed by:

